

FLIMSY PCR TUBES? – TRY ULTRAFLUX RIGID TUBES



Key Features:

- Angularly attached caps prevent cap lids and hinges from interfering with each other, allowing for easier operation in an SBS format rack.
- Robust interwell linkages make the tubes rigid, decreasing the chance of spillage.
- Unique two-material manufacturing creates the perfect environment for consistent, reliable qPCR performance. White wells offer enhanced fluorescence reflection and optically clear caps ensure that the fluorescent signal is passed through unaffected.
- Certified free of detectable RNase DNase, DNA, PCR inhibitors, and tested pyrogen-free.

IN STOCK NOW

| DESCRIPTION | CODE | PACK | PRICE |
|------------------------------------------------------------------------------|-------------|--------|----------|
| 0.2mL PCR Tubes UltraFlux I 8-Strip w Individually-Attached Flat Caps, Clear | SSIB3247-00 | 120pcs | \$106.75 |

Compatibility:

Applied Biosystems® 2700, 2720, 6100, 5700, PRISM 7000, 7300, 7500, 9600, GeneAmp® 9700, PRISM 310

Bio-Rad® / MJ Research® Gene Cyclor, iCycler, iQ5, Mini Gradient, MyCycler, MyiQ, PTC-100, PTC-200, PTC-225 Tetra, Dyad

Eppendorf® Mastercycler, Mastercycler ep Gradient, Mastercycler ep Realplex, Mastercycler Nexus /Gradient, Mastercycler Nexus Eco

Biometra® T1 Thermal Cyclor, Tgradient, Trobot, Uno, Uno II

| | | | |
|--------------------------------------------------------------------------------------------------------------------|-------------|--------|----------|
| 0.1 mL Low Profile, Strip of 8. PCR Tubes UltraFlux I 8-Strip Low-Profile w Individually-Attached Flat Caps, Clear | SSIB3248-00 | 120pcs | \$109.80 |
|--------------------------------------------------------------------------------------------------------------------|-------------|--------|----------|

Compatibility:

Applied Biosystems® 7500 Fast

Bio-Rad® / MJ Research® Chromo-4, Opticon / Opticon2

Eppendorf® Mastercycler, Mastercycler Nexus /Gradient, Mastercycler Nexus Eco

Promotion is valid until 30th August 2020 or while stocks last • Prices exclude GST and local freight

OTHER PCR TUBE OPTIONS:

IN STOCK NOW

| DESCRIPTION | CODE | PACK | PRICE |
|--------------------------------------------------------------------------------------------------------------------------------------------|-------------|---------|----------|
| PCR Tubes SnapStrip® 8-Strip w Individually- Attached Dome Caps, Clear. Single interwell linkage allows easy split in to smaller sections. | SSIB3240-00 | 120pcs | \$106.75 |
| PCR Tube 0.2ml Individual PCR Tube w Optically-Clear Flat Cap, Clear | SSIB3225-00 | 1000pcs | \$48.90 |
| PCR Tube 0.2ml Individual w Dome Cap, Clear | SSIB3210-00 | 1000pcs | \$60.50 |

Promotion is valid until 30th August 2020 or while stocks last • Prices exclude GST and local freight



Advantages of UltraFlux and SnapStrip Tubes:

Angular Cap Attachment


- Angularly attached caps on UltraFlux® i and SnapStrip® tubes prevent cap lids and hinges from interfering with each other, allowing for easier operation in an SBS format rack.

Ultra-Thin Walls

- Ultra-thin walls ensure efficient thermal transfer and maximum yield, yet are strong enough to prevent crushing.

Tube Color and Clarity

- Clear tubes provide high clarity for visualizing contents.
- White tubes are available with optically-clear caps, thanks to proprietary two-material manufacturing. The white wells offer enhanced fluorescence reflection and the optically-clear caps ensure that the fluorescence signal is passed through unaffected, creating the perfect environment for consistent and reliable qPCR performance.

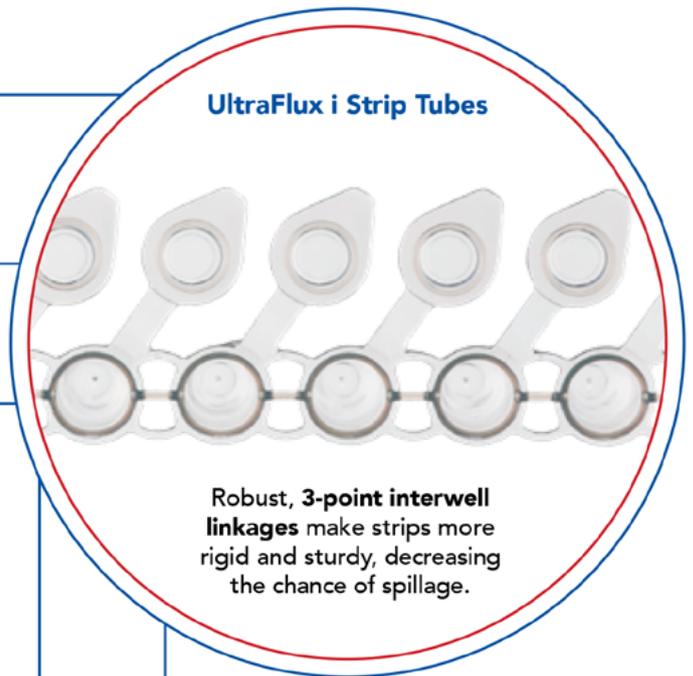


Optically-clear caps ensure that the fluorescence signal is passed through uninhibited.



SnapStrip single interwell linkage makes strips easy to split into smaller sections.

White tubes offer enhanced fluorescence reflection.

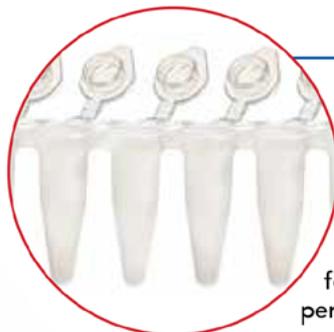

Unique two-material manufacturing creates the perfect environment for consistent, reliable qPCR performance.