


## JumboMix® 3500 WarmMix® CC

3500 mL lab blender

Ref 033 230



standards:

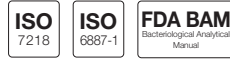

Manufactured under:



Product made for INTERSCIENCE by Interlab, an ISO 9001 certified company

**JumboMix® 3500 WarmMix® CC** enables the quick homogenizing of large samples. Avoiding all risks of cross-contamination, it is an easy-to-use and powerful lab blender. Its regulated door (from 20°C to 50°C) provides perfect conditions for the Trichinella testing in pork.

### TECHNICAL SPECIFICATIONS

- All stainless steel 304L
- Regulated heated door from 20°C to 50°C (68°F to 122°F)
- Multifunction digital screen
- Optimal blending volume: 200 - 3750 mL
- Variable speed (1.5 - 3 - 4.5 - 6 strokes/s)
- Variable blending time (1 s - 1 h or ∞ / count-down)
- **Adjustable blending power (patented)**
- **Side-by-side paddle stop**
- **Click & Clean® removable paddles (patented)**
- **Q-Tight® closing: tight closing for no more leaks from the top of the bag**
- Fully opening door (270°)
- Security drip tray
  
- Voltage and frequency: 100-127V~ / 220-240V~ 50-60Hz
- Fuses: T 4A L 240V / T 2A L 240V 5x20 mm
- Max power: 400 W
  
- Dimensions (w x d x h): 63 x 50 x 47 cm, weight: 63 kg
- Box (w x d x h) : 62 x 58 x 58 cm, weight: 61 kg
  
- 3-year warranty
- **Shock absorbers lifetime warranty**
  
- In compliance with: ISO 7218, ISO 6887-1 and FDA BAM (Bacteriological Analytical Manual)
- Manufactured under CE, RoHS, WEEE - Products made for INTERSCIENCE by Interlab, an ISO 9001 certified company
- Designed and made in France

### DELIVERED WITH

- 1 security drip tray
- 1 power cable
- 1 user's manual
- 1 quick user guide
- 1 free pack of blender bags

### OPTIONS / ACCESSORIES

- Ref 037 000 **JumboStand:** Ergonomic table for **JumboMix® 3500**